



Comfort through flexibility

Highly-adaptable System One humidification plus a Heated Tube option provide more ways for patients to get comfortable.

Building on Philips' highly successful System One humidification, Philips Respironics introduces an additional humidification solution that enhances performance, flexibility and comfort. With the System One Heated Tube, we now have a solution for exceptional patient comfort, even in cold environments.

There's no heated debate

System One with Heated Tube provides patients with;

- Air temperature control
- Improved humidification
- Rainout protection

There's little debating, the System One Heated Tube, combined with proven System One humidification, is a win-win scenario for patients and providers alike.

PHILIPS

RESPIRONICS

sense **and** simplicity

Get warm and cozy with System One Heated Tube

It's just the right combination of warmth and humidity

The System One device incorporates a room temperature and humidity sensor with a temperature sensor at the end of the Heated Tube to deliver the selected level of humidity to the patient while providing protection against rainout. As a result, higher levels of humidity can be achieved for patients who could benefit from humidification.

With System One Heated Tube, you get increased humidification plus the rainout protection you've come to trust and expect from Philips Respironics.




A top-notch design

System One Heated Tube was designed with:

- top-mounted connection for easy access
- durable, simple one-click connection
- 270-degree swivel mount

More than a number

Humidification is now more personalized with System One Heated Tube. It allows patients to select their preferred humidity and temperature levels throughout the night, regardless of changing room conditions.

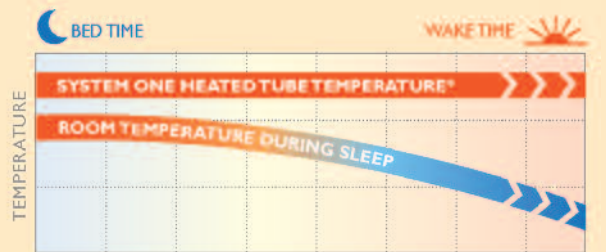

* 5 temperature settings available

Ordering information	Part number
Heated humidifier standard	IN6H
Heated tube humidifier	IN6T
Upgrade kit	6T-CONVKIT
System One Heated Tube	SYS1HT15
Water chamber	1063785
80 watt power supply	1091399
Power cord Europe	1039014
Power cord UK	1039026

Respironics, REMstar, System One, and A-Flex are trademarks of Koninklijke Philips Electronics N.V. and its affiliates. All rights reserved.

Please visit www.philips.com/systemone



© 2012 Koninklijke Philips Electronics N.V. All rights are reserved.

Philips Healthcare reserves the right to make changes in specifications and/or to discontinue any product at any time without notice or obligation and will not be liable for any consequences resulting from the use of this publication.

CAUTION: US federal law restricts these devices to sale by or on the order of a physician.

Broudy NA XX/XX/XX MCI 4104926 PN 1101231

Philips Healthcare is part of Royal Philips Electronics

Philips Respironics International Headquarters
+33 1 47 28 30 82
Philips Respironics Asia Pacific
+65 6882 5282
Philips Respironics Australia
+61 (2) 9947 0440
Philips Respironics United Kingdom
+44 800 1300 845

Philips Respironics
1010 Murry Ridge Lane
Murrysville, PA 15668

Customer Service
+1 724 387 4000
800 345 6443 (toll free, US only)
www.philips.com/respironics